

Specification Sheet

FibroLife™ S2 Cell Culture Medium

LL-0011

LM-0001

LS-1038



This product is for Research Use Only.
This product is not approved for human or
veterinary use or for use in in vitro
diagnostics or clinical procedures.



LIFELINE
CELL TECHNOLOGY

Better Solutions for
Breakthrough Results

FibroLife S2 Cell Culture Medium

FibroLife S2 is a variation of our FibroLife cell culture medium optimized for low serum culture of Fibroblasts and Mesenchymal Stem Cells-Pre-Adipocytes. FibroLife S2 supports the growth of these cells in a 2% serum environment at rates equal to or greater than comparable media supplemented with 2 to 10% FBS.

FibroLife Basal Medium contains no antimicrobials and no phenol red. Lifeline offers antimicrobials and phenol red; however, they are not required for eukaryotic cell proliferation.

FibroLife S2 Cell Culture Medium Components

FibroLife S2 is offered in a kit format composed of basal medium and associated supplements and growth factors called “LifeFactors™.” This kit allows you to prepare fresh medium in your laboratory, extending shelf life and enhancing performance. FibroLife Basal Medium is provided in a light-protected 500 mL bottle. The remaining LifeFactors (growth factors and supplements) are packaged in a convenient gas-impermeable pouch for easy storage.

Product	Part No.	Volume	Final Concentrations in Supplemented Medium	Storage
FibroLife S2 Medium Complete Kit, (FibroLife Basal Medium, FibroLife S2 LifeFactors Kit)	LL-0011			2-8°C when prepared
FibroLife Basal Medium	LM-0001	480 mL		2-8°C
FibroLife S2 LifeFactors Kit	LS-1038			-20°C
rh FGF basic LifeFactor	LS-1002	0.5 mL	5 ng/mL	-20°C
rh Insulin LifeFactor	LS-1004	0.5 mL	5 µg/mL	-20°C
Ascorbic Acid LifeFactor	LS-1005	0.5 mL	50 µg/mL	-20°C
L-Glutamine LifeFactor	LS-1006	18.75 mL	7.5 mM	-20°C
Hydrocortisone Hemisuccinate LifeFactor	LS-1007	0.5 mL	1 µg/mL	-20°C
FBS LifeFactor	LS-1012	10 mL	2%	-20°C
Antimicrobial Supplement: Gentamicin and Amphotericin B (Provided with purchase of LL-0011)	LS-1104	0.5 mL	Gentamicin 30 µg/mL Amphotericin B 15 ng/mL	-20°C
Optional Supplements	Part No.	Volume	Concentrations of Supplement	Storage
Phenol Red Supplement	LS-1009	1 mL	33 mM	RT

Quality Testing for Guaranteed Consistency and Reproducible Results

Lifeline Cell Technology manufactures products using the highest quality raw materials and incorporates ISO style quality assurance in every production run. Exacting standards and production procedures ensure lot-to-lot consistency. Every production lot of FibroLife S2 Medium is extensively tested using human fibroblasts. FibroLife S2 Medium and Gentamicin-Amphotericin B ([LS-1104](#)) are cell culture tested separately.




Call Lifeline Technical Service and Sales at **877.845.7787**

Or visit lifelinecelltech.com for more information

Lifeline Cell Technology | 8415 Progress Drive, Suite T | Frederick, MD 21701

©2024 Lifeline Cell Technology, an International Stem Cell Company | All Rights Reserved



- Sterility Testing Negative for bacterial and fungal growth
- pH 7.4 +/- 0.2
- Cell Testing Rate of proliferation and morphology
- Osmolality 270 +/- 10 mOsm
- Endotoxin < 0.5 EU/mL

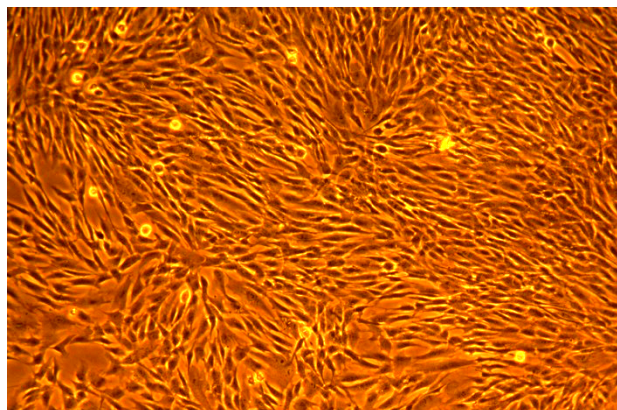

FibroLife™ S2 medium grows Fibroblasts and HMSC-Pre-Adipocytes in low serum conditions at rates that meet or exceed other commercially available serum-containing media while maintaining excellent cell morphology. In comparisons with other commercially available media, FibroLife S2 medium shows superior proliferation at different seeding densities.

Shown left: Human Dermal Fibroblasts-Neonatal, passage 2, 5 days after inoculation with 2,000 cells/cm² (100X).

Shown below: Human Dermal Fibroblasts-Neonatal inoculated at low seeding densities in 24-well plates. Higher fluorescence represents a greater number of cells per well and therefore better proliferation.

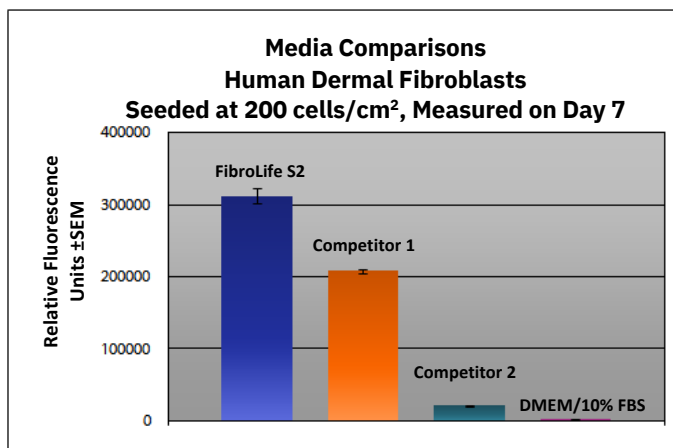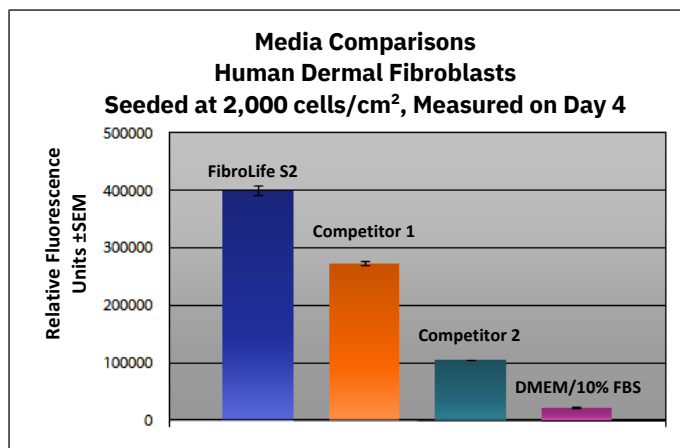

The Lifeline® Guarantee

Lifeline's rigorous quality control ensures sterility and performance to standardized testing criteria. Upon request, Lifeline will provide lot-specific QC test results, material safety data sheets and certificates of analysis. See complete guarantee/warranty statement at lifelinecelltech.com or contact your Lifeline representative for more information.



Call Lifeline Technical Service and Sales at **877.845.7787**

Or visit **lifelinecelltech.com** for more information

Lifeline Cell Technology | 8415 Progress Drive, Suite T | Frederick, MD 21701

©2024 Lifeline Cell Technology, an International Stem Cell Company | All Rights Reserved